

karma



Ergo Recline


Reclining Transfer Wheelchair

The Ergo Recline offers a near-horizontal recline with flip-back detachable armrests for convenient lateral transfers and excellent comfort.

Patented Multi-functional Armrests

Flip-back, height-adjustable, and fully detachable armrests accommodate a wide range of use cases, such as changing clothes and transferring.



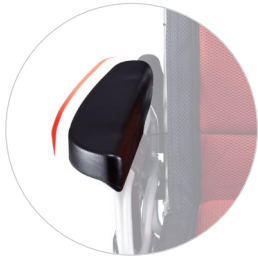

Boat Armpad

Wider than traditional armpads and featuring a comfortable, contoured shape, they are easily re-positioned to accommodate self-propulsion and transport.



Ergonomic Headrest with 3D Pillow

Provides lateral and occipital support to relieve neck and shoulder pressure without restricting vision. Velcro backing allows you to adjust the headrest position to create the perfect fit.



Elevating Legrests

Angle-adjustable, with calf pads, for improved comfort and support.

Hydraulic Pressure Support System

The Ergo Recline's recline system contains two pneumatic gas struts, which provide a smooth, controlled tilting motion and actively assist the attendant return the backrest to the upright position.

Ergo Recline Specifications

Seat Width* Seat Depth (inch)	16*16 / 18*16
F/R wheels (inch)	8 / 14 、 8 / 20
Seat Surface Height at Front Edge (mm)	500
Backrest Height (mm)	495
Armrest Height (mm)	220 - 300
Calf Length (mm)	365 - 440
Seat to Backrest Angle (X°)	8 - 80
Max. User Weight (kg)	115
Total Mass (kg)	21.1 / 21.2 / 21.8 / 21.9
Overall Length with legrest (mm)	1190 - 1815
Overall Width (mm)	590 / 640 / 690
Overall Height (mm)	700 - 1300
Folded Length (mm)	790 / 850
Folded Width (mm)	340 / 400
Folded Height (mm)	940

Remark : Karma reserves the right to modify information herein without further notice.
 * / " means optional ; " ° " means adjustable ±1cm ±0.5kg
 Rev.00 Jul. 2024



With a 161° reclining angle, the seating pressure can be redistributed across the entire body.



Scan to learn more about Ergo Recline



KARMA MEDICAL PRODUCTS CO., LTD.
 T +886-5-206 6688
 F +886-5-206 7788
 Email : globalsales@karma.com.tw