



Affordable, Durable, Stylish



Soma Agile

Configurable, with standard flip-up armrests and swing away footrests, to meet a wide range of needs.



Soma 200 Series

Ultra lightweight for travel with patented swing away footrests for easy transfer.



Soma 100 Series

Economic aluminum wheelchair, featuring fixed armrests and footrests.

soma

Designed and Made by Karma

SOMA Agile



Designed with both users and attendants in mind, Agile is the new benchmark for standard wheelchairs. Detachable, swing-away footrests and flip-up armrests make transfers quick and easy, while the aircraft-grade aluminum frame is lightweight and compact to maximize portability and ease of use.

A wide range of options and accessories allow Agile to fulfill the needs of a wide range of users.



Rear Wheel Extension Plate (Optional)

Adjust the horizontal position of the rear wheel to adjust propelling efficiency or stability, according to user needs.



Tension Adjustable Backrest (Optional)

The internal structure of the backrest is formed by a series of individually-adjustable Velcro straps, which can create a user-specific fit to maximize comfort and support.



Quick-Release Rear Wheels

The rear wheels are removable with the push of a button, minimizing Agile's size for storage and transport.



| Rear Wheel | Seat Dimensions (W x D) | |
|-----------------|-------------------------|-----------|
| Self-propel 24" | 16" x 16" | 16" x 18" |
| | 18" x 16" | 18" x 18" |
| | | 20" x 18" |
| Self-propel 22" | 16" x 16" | 16" x 18" |
| | 18" x 16" | 18" x 18" |
| Transport 14" | 16" x 16" | 16" x 18" |
| | 18" x 16" | 18" x 18" |
| | | 20" x 18" |

FRAME COLOR





Easy Transfer with Agile

Standard flip-up armrest and swing away, detachable footrests make transferring in and out of the wheelchair a breeze.

Optional upgrades include, height-adjustable armrest, elevating legrests, quick release wheels, and more!




One Frame, Three Rear-Wheel Sizes

Agile's frame is compatible with 14", 22", and 24" rear wheels; converting a transport chair to a self-propelled wheelchair is as simple as adjusting the conversion plates and switch the rear wheel.



SOMA 200 Series

Lightweight aluminum wheelchair with patented swing away footrests; at only 10.2kg, the SOMA 215 is perfect for daily transport.

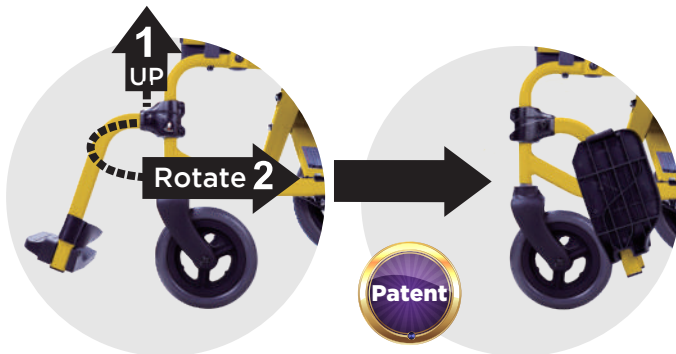

Never Lose a Footplate Again!

The patent swing away footrests are easy to swing away without the need of detaching it, making everyday transfers, storage, and transportation smooth and worry-free.



Seamless Elliptical Frame

Crafted using minimal components and welds, resulting in both a refined aesthetic and exceptional durability.



Detachable Upholstery

Intuitive design makes it easy to detach and clean the upholstery.



Ergonomic Push Handles

Offers a stable and comfortable operating experience for attendants.



Attendant Tipping Lever

Allow caregivers to overcome obstacles with ease.


FRAME COLOR



Bright Silver White



Yellow

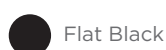

Flat Black

Various wheel sizes and colors to suit your needs and style!



Low Seat to Floor Height

The accessible 43 cm seat to floor height enables most users to propel the wheelchair with their feet on the floor. (SOMA 215TL 48 cm)



Self-propel

| | Caster x Rear wheel size | Seat to Floor Height | Seat Dimensions (W x D) |
|------------|--------------------------|----------------------|-------------------------|
| SOMA 215 | 6" x 20" F | 43 | 16" x 15" 17" x 15" |
| | 6" x 14" F | 43 | 16" x 15" 17" x 15" |
| SOMA 215TL | 7" x 14" F | 48 | 18" x 17" |
| | 7" x 20" F | 48 | 18" x 17" |
| | 7" x 22" Q | 48 | 18" x 17" |



Transport



Anti-Slip Footplates

Improved grip allows users to better maintain proper foot position.

| | SOMA 105 | SOMA 215 | SOMA 215TL | Agile |
|----------------------|----------|----------|------------|-----------|
| Foldable Backrest | ✓ | ✓ | ✓ | ✓ |
| Swing-away Footrests | | ✓ | ✓ | ✓ |
| Flip-back Armrests | | | | ✓ |
| Anti-tippers | | ✓ Option | ✓ Option | ✓ |
| User Height (cm) | 162cm ↑ | 162cm ↓ | 162cm ↑ | 162cm ↑ |
| Total Weight (kg) | 11kg | 10.2kg | 10.9kg | 12-14.1kg |

SOMA 105

Lightweight aluminum wheelchair featuring fixed full-length armrests and fixed footrests with flip-up footplates



| Rear Wheel | Seat Dimensions (W x D) |
|-----------------|-------------------------|
| Self-propel 22" | 16" x 17" 18" x 17" |
| Transport 16" | 16" x 17" 18" x 17" |

Full-Length Armrest Design

Offers ample support for transfers and repositioning.



16" transport rear wheel option



SOMA 100 Series Features



Ergonomic Push Handles

Offers a stable and comfortable operating experience for attendants.



Attendant Tipping Lever

Allow caregivers to overcome obstacles with ease.



Foldable Backrest

For easy transport and storage.

SOMA 103

Lightweight aluminum wheelchair featuring fixed desk-length armrests and fixed footrests with flip-up footplates



| Rear Wheel | Seat Dimcnsions (W x D) |
|-----------------|-------------------------|
| Self-propel 22" | 17" x 17" |
| Transport 16" | 17" x 17" |

Desk-Length Armrests

Allow the wheelchair user to move closer to desks and tables.



16" transport rear wheel option



Detachable Upholstery

Intuitive design makes it easy to detach and clean the upholstery.



Anti-Slip Footplates

Improves foot positioning and prevents slipping.

| Product Name | Nominal seat width (mm)* | Nominal seat depth (mm)** | Front caster / Rear wheel size (in) | Seat surface height at front edge (mm) | Backrest height (mm) | Armrest to seat distance (mm) | Footrest to seat distance (mm) | Mass of the heaviest part (kg) | Total mass (kg) |
|--------------|--------------------------|---------------------------|-------------------------------------|--|----------------------|-------------------------------|--------------------------------|---|---|
| SOMA 215 | 425 | 395 | 6 x 14 | 430 | 410 | 220 | 390 - 440 | 10.2 | 10.2 |
| | 425 | 395 | 6 x 20 | 430 | 410 | 220 | 390 - 440 | 11.2 | 11.2 |
| SOMA 215TL | 460 | 430 | 7 x 14 | 480 | 430 | 220 | 430 - 474 | 10.5 | 10.9 |
| | | | 7 x 20 | | | | | 11.7 | 12 |
| | | | 7 x Q22 | | | | | 8 | 12.4 |
| SOMA 100 | 400 | 430 | 7 x 16 | 500 | 390 | 220 | 460 - 510 | 10.2 | 10.4 |
| | 450 | | | | | | | 10.5 | 10.7 |
| | 400 | 430 | 7 x 22 | 480 | 390 | | 395 - 445 | 11.1 | 11.3 |
| | 450 | | | | | | | 11.4 | 11.6 |
| SOMA 105 | 400 | 430 | 7 x 16 | 500 | 390 | 220 | 460 - 510 | 10.3 | 10.5 |
| | 450 | | | | | | | 10.6 | 10.8 |
| | 400 | 430 | 7 x 22 | 480 | 390 | | 395 - 445 | 11.2 | 11.4 |
| | 450 | | | | | | | 11.5 | 11.7 |
| SOMA 103 | 420 | 430 | 7 x 16 | 500 | 400 | 235 | 395 - 420 | 10.4 | 10.4 |
| | | 430 | 7 x 22 | 500 | 400 | 235 | 395 - 420 | 11.6 | 11.6 |
| AGILE | 16 | 16 | 7 x 14 | 480 | 430 / 510 | 220 | 380 - 480 | 11.1 - 11.6 11.4 - 11.9 11.6 - 12.1 | 12.3 - 12.8 12.6 - 13.1 12.8 - 13.3 |
| | 18 | 16 | | 480 | | 220 | | | |
| | 16 | 18 | | 480 | | 220 | | | |
| | 18 | 18 | | 480 | | 220 | | | |
| | 20 | 18 | | 480 | | 220 | | | |
| | 16 | 16 | 7 x 22 | 480 | 430 / 510 | 220 | 380 - 480 | 11.7 - 12.2 12.0 - 12.5 12.2 - 12.7 | 12.9 - 13.4 13.2 - 13.7 13.4 - 13.9 |
| | 18 | 16 | | 480 | | 220 | | | |
| | 16 | 18 | | 480 | | 220 | | | |
| | 18 | 18 | | 480 | | 220 | | | |
| | 20 | 18 | | 480 | | 220 | | | |
| | 16 | 16 | 8 x 14 | 505 | 430 / 510 | 220 | 380-480 | 11.0 - 11.5 11.3 - 11.8 11.5 - 12.0 | 12.2 - 12.7 12.5 - 13.0 12.7 - 13.2 |
| | 18 | 16 | | 505 | | 220 | | | |
| | 16 | 18 | | 505 | | 220 | | | |
| | 18 | 18 | | 505 | | 220 | | | |
| | 20 | 18 | | 505 | | 220 | | | |
| | 16 | 16 | 8 x Q24 | 505 | 430 / 510 | 220 | 380-480 | 8.6 - 9.1 8.9 - 9.4 9.1 - 9.6 | 13.9 - 14.4 14.2 - 14.7 14.4 - 14.9 |
| | 18 | 16 | | 505 | | 220 | | | |
| | 16 | 18 | | 505 | | 220 | | | |
| | 18 | 18 | | 505 | | 220 | | | |
| | 20 | 18 | | 505 | | 220 | | | |

The actual product specifications differ from the data in the table by ±1 cm ±0.5 kg.



A Better Wheelchair, A Better Fit

KARMA MEDICAL PRODUCTS CO., LTD.
T +886 5 206 6688 F +886 5 206 7788

globalsales@karma.com.tw
www.karmamedical.com



*Karma reserves the right to modify information herein without further notice. • Rev. 01 Aug. 2023