



FDA MDA  CE ISO 9001

104-40585 Rev. 02 Oct. 2014

Discover the path and make it yours.

**karma**  
mobility solutions



No. 2363 Sec. 2 Da-Shiue Rd.,  
Min-Hsiung Shiang, Chia-Yi 621, Taiwan  
Tel: 886-5-206-6688 Fax: 886-5-206-7788  
E-mail: [globalsales@karma.com.tw](mailto:globalsales@karma.com.tw)

[www.KarmaMedical.com](http://www.KarmaMedical.com)

**karma**  
mobility solutions

KARMA MEDICAL PRODUCTS CO., LTD.

**MVP502**

Most Value Performance *KM-5000.2*




# MVP 502

Most Value Performance

## After reclining, what do you not feel?

No matter you sit or recline, MVP 502 is contoured to fit your body, relieving the pressure and shear force.

The revolutionary coupling system of MVP 502 offers both the user and attendant an ultimate comfort while reclining. Moreover, the innovated anti-sliding design and the V-shaped seating can effectively distribute pressure for relief of bottom, greatly reducing the occurrence of back sores, and keeping seating position really pleasurable.

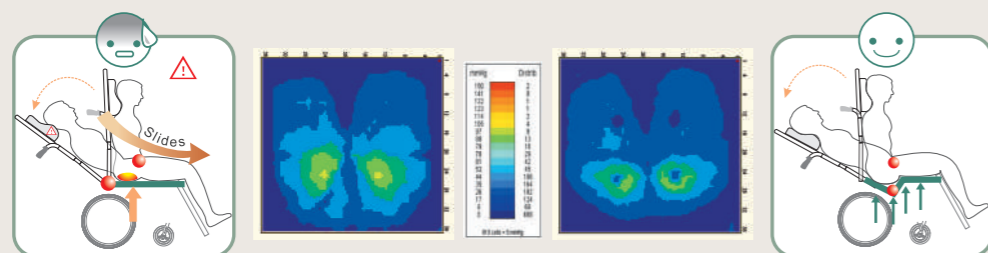

## New Patented Anti-sliding V-seating system

The V-shaped seat smoothly shifts the bottom 20° downward, minimizing shear force when the backrest reclines 160°. By doing this, it not only prevents sliding, but also effectively relieves pressure on the bottom, which reduces the risk of decubitus ulcers.

## Features




### Rigidifying Push-bar (Lateral Support bar)

The rigidifying push-bar strengthens back support, offering greater rigidity to the chair, and making it easy to push and steer the chair. Folding structure for added convenience.



### Multi-functional ergonomic armrest

The backrest-armrest coupling structural design allows the armrest to switch position ergonomically for user's elbow when the backrest reclines from 93° to 160°. It is also height-adjustable and easy to flip back.

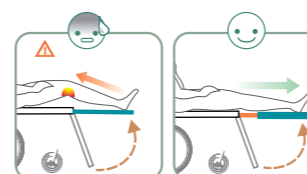

### Axle-traveling footrest

The elevating swing-away footrest automatically stretches during elevation, so that the user's legs can extend smoothly without compression, which effectively eases the pressure concentration at abdomen and offers more extending space for users with longer legs.



### Strap-tension adjustable backrest and seat

Unlike tradition, the back- and seat-strap design provides unrestricted adjustments according to user's condition while stabilizing the sitting position, preventing spinal deformity and kyphosis.



### 3D elevating legrest support

Ergonomically designed, allowing greater micro-adjustment: up/down, left/right, inward/ outward.



### Elite gas spring

The high performance hydraulic gas spring aids upward push during sitting up and reduces shock during reclining; it achieves dual shock reduction by enhancing the stabilization of comfort reclining.



### Ergonomic headrest

Improving head position and enhancing greater relaxation, the ergonomically designed headrest can be adjusted to relieve pressure for shoulder and neck, and offers better vision for attendants.



### Carefree pressure-relieving seat upholstery

The detachable backrest/seat upholstery is 3 cm thick, which provides a fluffy sensation offering both great comfort and ventilation.



### Foldable and Compact

It folds easily and compact after folding, which makes it easy for storing and traveling.



AEGIS Microbe Shield permanently bonded anti-microbial barrier upholstery protects from odor, staining and deterioration from bacteria, fungus and other microorganisms. It is a shield for your health.

## Optional accessories



Armrest, ergonomic model, offers greater comfort and fitness for elbow (for 14" rear wheel only)



Quick release 22" · 20" rear wheels

■ Multiple seat width selection via customization: 15" X 16", 20" X 18"